

RELAX

FORM FIT BACK



Rear Back Angle



Removable Cube Foam

The **Relax Form Fit** is an aluminum light weight formed back with hardware that allows for multiple adjustments. The Relax Form Fit with its wide range of laterals is suited for complex seating and positioning needs. Cubes are removable to custom fit for spinal anomalies.

Quick Reference Guide

RELAX SERIES - FORM FIT BACK

REGULAR SIZES ADP: [SEBWM0026](#)

Width = 16" - 18"

Part Number: Refer below

Height = 18"

CUSTOM ADP: WIDTH [SEBWM0035](#), HEIGHT [SEBWM0040](#)

Width = 15" or <, 19" or >

Part Number: Refer below

Height = 17" or <, 19" or >

"<" means less and ">" means more

PART NUMBERS

SHUF (W X H) **No Laterals**

SHUS (W X H) **Gentle Laterals**

SHUM (W X H) **Snug Laterals**

SHUD (W X H) **Full Laterals**

Fill in width (W) and height (H) to get part number.

E.g. For a W = 15" and H = 17" Back; Part number would be SHU1517

Poly Foam Cubes Part Numbers:

QUB2x2


Poly Foam Cubes 2x2

QUB3x2

Poly Foam Cubes 3x2


RELAX

RELAX FORM FIT COMES WITH THE FOLLOWING BACK SHELL OPTIONS

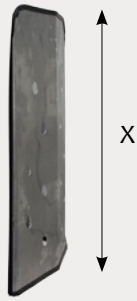


NO LATERAL DIMENSIONS

X = Standard shell height 18"

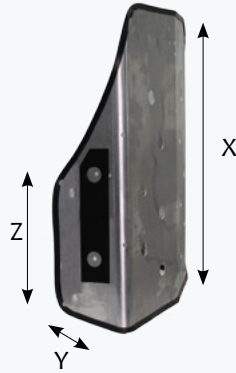


GENTLE SHELL LATERAL DIMENSIONS

X = Standard shell height 18"
Y = Lateral depth 3" finished
Z = Lateral height 10"

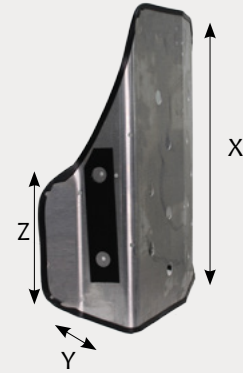

FULL SHELL LATERAL DIMENSIONS

X = Standard shell height 18"
Y = Lateral depth 5" finished
Z = Lateral height 10"

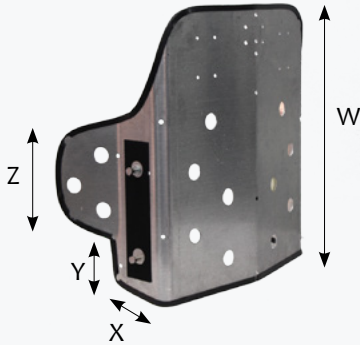


SNUG SHELL LATERAL DIMENSIONS

W = Standard shell height 18"
X = Lateral depth 5" finished
Y = Space below lateral 4"
Z = Lateral height 6" above cutout

KEY FEATURES

- Wide range of laterals with removable cubes for custom fit
- Comfortable High Resilience Foam
- 3 lbs. Aluminum, Light Weight Shell
- Pre Punched holes to mount headrest and laterals
- Adjustable/Removable Lumbar Support
- Easy to Install and Remove
- Height Adjustable

OVERLAYS

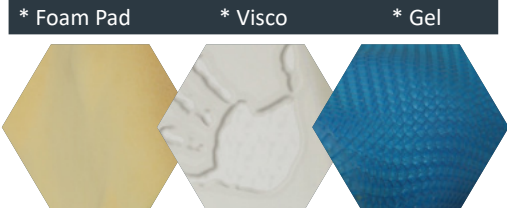

ADP: SEBND0015 SEBND0015 SEBND0020

* Up-charge items

SPECIFICATIONS	RELAX SERIES - FORM FIT BACK
ADP	SEBWM0026, CUSTOM W SEBWM0035, H SEBWM0040
STANDARD WIDTH RANGE	16" - 18"
CUSTOM WIDTH RANGE	15" or less and 19" or more
STANDARD HEIGHT RANGE	18"
CUSTOM HEIGHT RANGE	17" or less and 19" or more
AVAILABLE HARDWARE STYLES	One style which is universal to all our backs
ANGLE ADJUSTABLE	25° anterior adjustability, 20° of posterior adjustability
WIDTH ADJUSTABLE	2" of adjustability through the pins on the back
FITS TUBING SIZE	7/8" tubing or 1"
HEADREST MOUNT	Backs pre-drilled for headrest supports and lateral supports