

# NV TILT

TILT WHEELCHAIR




ADP CODE: WM0001076


The **NV TILT** is a light weight tilt chair designed for adjustment and maneuverability. The **NV TILT** tilt system provides a change of orientation to redistribute pressure from one area to another.



Center Pivot Swing In And Swing Out  
Front Rigging  
(Optional)

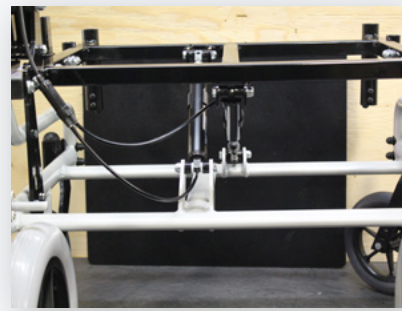

Dual Cylinders Allows For Consistent  
Tilt For All Weight Capacities

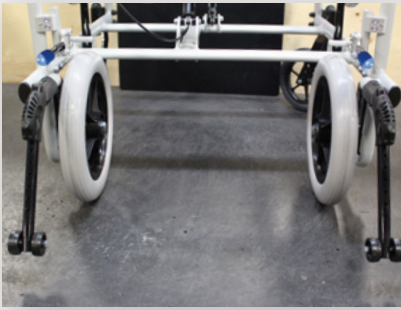

Narrow Configuration  
(Inside Wheel Mount)



Rear Axle Adjustability



Rear Depth Adjustable Frame

## ABOUT NV TILT

The **NV TILT** is designed for adjustment and maneuverability, with comfort being another major point in this design. Our **NV TILT** will be your choice in tilt wheelchairs and we guarantee this product will exceed your expectations. The **NV TILT** offers a wide range of dynamics that allow this wheelchair to accommodate many different situations, assuring a right fit for many, even in the most difficult positions.

## CUSTOMIZE ANYTHING

Power Plus Mobility's advanced machinery specializes in manufacturing custom parts to customers' specifications from either a drawing or a sample. The parts we manufacture can be simple or more complex and can be manufactured in small or large quantities. Contact your local sales representative with your specifications.

# NV TILT



## STANDARD FEATURES

- Light Weight Frame
- Depth Adjustable Frame
- ABS Seat Pan
- Dual Gas Struts For Easy Tilt (Up To 55° )
- Tilt Locked out at 33°
- New Recessed Lower Frame For Foot Propulsion
- Dual Gas Cylinders for Smooth & Easy Tilting
- Soft Set Gas Cylinder Release Mechanism
- Silver Grey or Black Frame Colour
- T-Style Arms
- Adjustable Stroller Handle
- Adjustable Caster Journals
- Recessed Axle Plates
- One Piece Elevating System
- Incontinent Seating
- Overall Wheelchair Width: 28" (18" Seat Width)

SPECIFICATIONS	NV TILT
SEAT WIDTH:	14" TO 36"
SEAT DEPTH:	16" TO 24"
SEAT HEIGHT:	(ADJUSTABLE): 12.25" TO 18.25"
TILT ANGLE:	0° TO 55°
BACK ANGLE:	(ADJUSTABLE): 0° , 5° , 10° OR 15°
ADJUSTABLE STROLLER HANDLE:	0° TO 180°
WEIGHT CAPACITY:	UP TO 600 lbs.
WEIGHT OF CHAIR:	60 lbs. (APPROXIMATELY)
	TWO-WAY LOCKING GAS STRUT
	SINGLE RELEASE, DUAL CYLINDERS
	STANDARD LOCKOUT COLLARS (33°)

NV TILT  
TILT WHEELCHAIR



Learn more about this product by scanning the QR Code